

# VACUUM INTEGRATED PUMP





## VIP



### DESCRIPTION

The ALPS™ VIP is an air expulsion integrated pump for elevated vacuum for prosthetic socket. The use of VIP provides patients continuous total contact when donning the prosthesis. This is essential for an even weight distribution and to distribute the suspension of the prosthesis over the whole surface of the residual limb. Elevated vacuum created by the VIP stabilizes the residual limb, overcoming all issues connected with the residual limb changing volume during the day.

### ADDED BENEFITS

-  Efficiently deals with stump volume fluctuation which in turn leads to improved prosthesis control and comfort
-  Reduces the intensity of peaks in pressure on the patient's limbs, improving comfort and providing greater freedom of movement
-  Helps preserve residual limb skin health by improving circulation and consequently increasing limb oxygenation
-  Does not affect the cosmetic appearance of the prosthesis




MAKING LIVES BETTER



#makinglivesbetter




# AVAILABLE IN 2 DIFFERENT VACUUM LEVELS

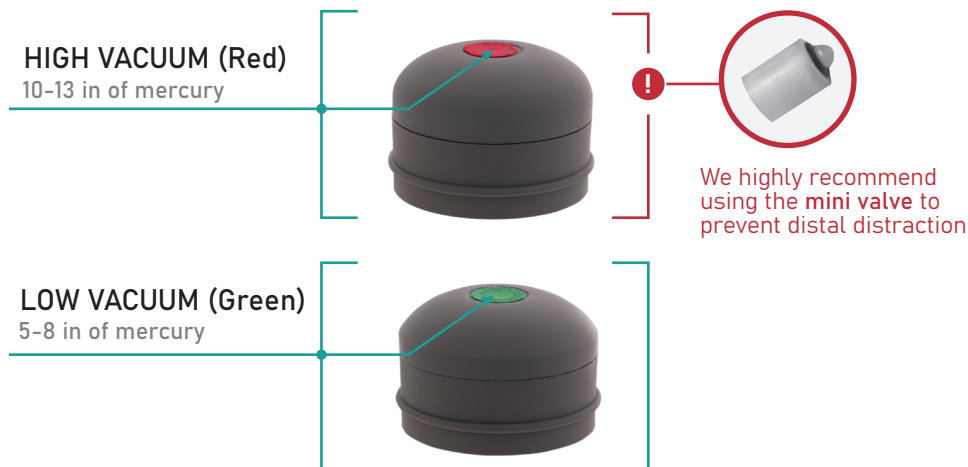


**CHOOSE LEVEL BASED ON:**

- ACTIVITY LEVEL
- SKIN CONDITION

## WITHOUT MOUNTING PLATE

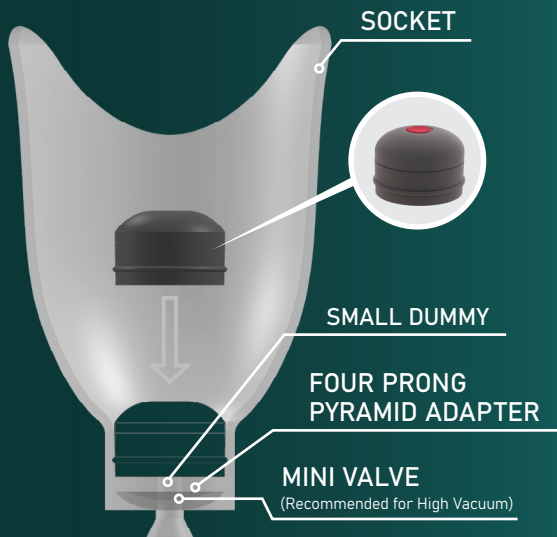

## WITH MOUNTING PLATE

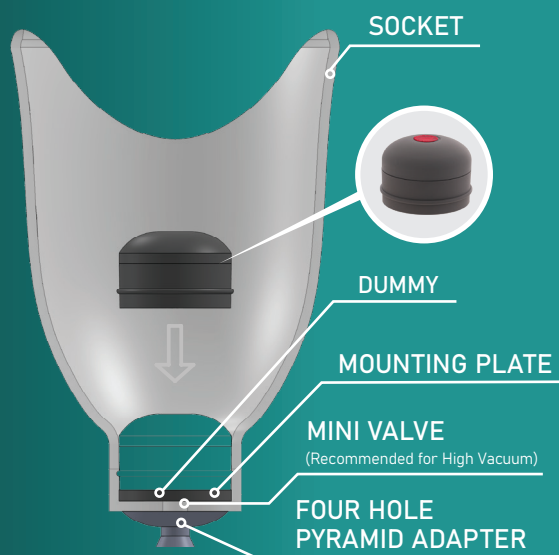

### PRODUCT SPECS

	HIGH VACUUM (RED)	LOW VACUUM (GREEN)
PART #	VIP - DH	VIP - DL
DIAMETER	66MM	66MM
HEIGHT	48MM	48MM
WEIGHT	75G	64G
MATERIAL	SILICONE	SILICONE

### USE CODES BELOW WHEN ORDERING

PART NAME	PART #
DUMMY	VE319 - D
SMALL DUMMY	VE319 - SD
MOUNTING PLATE	VE819 - P
LOCK (SMALL)	VIP-L1120-S
MINI VALVE (Highly Recommended)	MV-009
MINI VALVE KIT	LMV-009



Tel +420 377 223 127 | Fax +420 377 223 010

www.easyliner.eu | info@easyliner.eu

Božkovské náměstí 17/21, 32600 Plzeň - Czech Republic

