

# Designed for *My* Real Life



**Stephen | RV Instructor, CA**  
I'm on uneven terrain with holes when I camp, I used to worry about that, but with Kinnex I don't have to worry about that anymore. I can launch my boat **without fear** of stepping in the water. All the things that have prevented me from enjoying myself Kinnex removed them. **I love it.** With Kinnex I feel complete again.



**Mike | Retired Photographer, CA.**  
I love Kinnex because I can **get my life back** again, I can do the things I couldn't do like hiking and fishing and bike riding. **Life is better again.** With Kinnex I can get to places as a photographer I wasn't able to get to before, I'm able to go to those places with **confidence** and at this point I feel like I could do anything.



**Ramin | Property Manager, OK**  
Kinnex has helped me change my overall **quality of life.** You get the full motion, the **natural motion** that you miss with a fixed ankle. It takes a lot of the pressure and shock off of my actual joints. I feel **less fatigue** at the end of the day and my overall **quality and comfort** of my limb is tenfold better than it ever has been.



**Ellen | Rancher, CO.**  
Kinnex allows me to go on the uneven ground and I am **confident.** I don't worry about falling and I'm not having to pay much attention to what I'm doing. I can focus on the animals I'm taking care of and just live my life.

## Ordering made easy

### 1 Verify User's Measurement


- Patient weight
- Foot Size (24-30cm)
- Clearance requirements (see build height chart)
- Right or Left
- Sandal Toe/ Standard Toe

Foot Size	Build Height
24-25cm	7.2" / 184mm
26-28cm	7.3" / 186mm
29-30cm	7.3" / 186mm

### 2 Select Stiffness Category

Weight (lbs)	Impact Level		Weight (kg)
	Low	Moderate	
100-115	1	1	44-52
116-130	1	2	53-59
131-150	2	3	60-68
151-170	3	4	69-77
171-195	4	5	78-88
196-220	5	6	89-100
221-255	6	7	101-116
256-275	7		117-125

### 3 Complete Foot Module Part Number

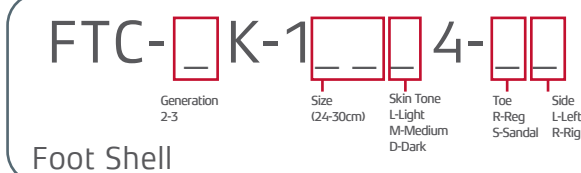

### Focused Solutions

Freedom Innovations is solely focused on developing world class lower limb solutions in close collaboration with prosthetists and people with amputations. These solutions encompass technological innovations, service, training, educational resources and consultative support to help ensure user satisfaction.


[www.freedom-innovations.com](http://www.freedom-innovations.com)    [www.freedomkinnex.com](http://www.freedomkinnex.com)

**Global Headquarters**  
3 Morgan  
Irvine, CA 92618  
toll free: 888-818-6777  
phone: 949-672-0032  
fax: 949-672-0084  
[www.freedom-innovations.com](http://www.freedom-innovations.com)

**European Head Office**  
Jaargetijdenweg 4  
7532 SX Enschede  
The Netherlands  
toll free: 00800 - 2806 2806  
phone: 31 (0) 53 20 30 300  
fax: 31 (0) 53 20 30 305  
[www.freedom-innovations.eu](http://www.freedom-innovations.eu)

**USA Manufacturing & Returns**  
425 East 400 North  
Gunnison, UT 84634

Connect with us  
  
 #freedomeveryday

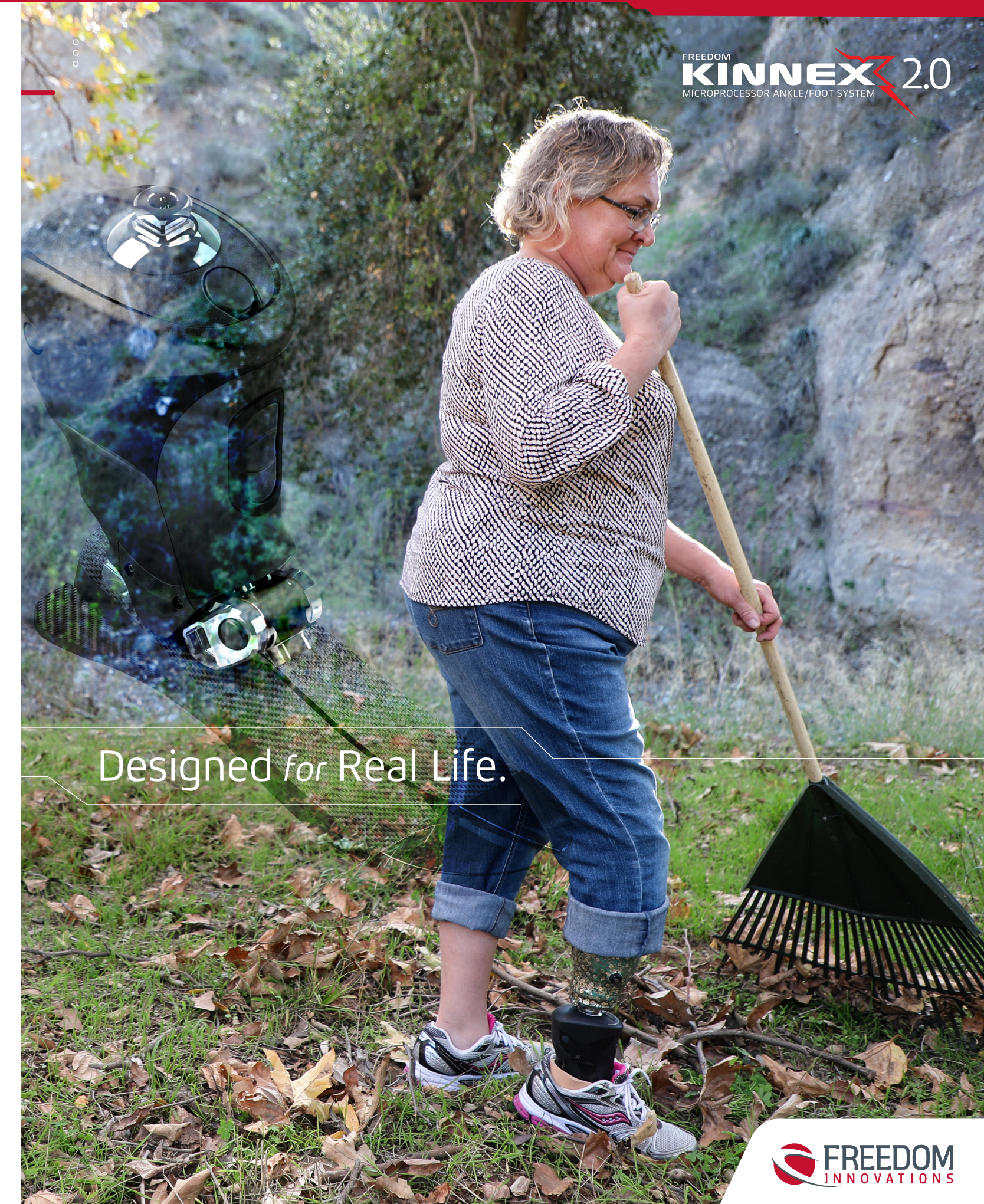

Designed for Real Life.



FREEDOM  
**KINNEX 2.0**



**Designed for real life.** Kinnex microprocessor ankle combines sophisticated microprocessor ankle innovation with Freedom's carbon fiber foot technology, providing users a uniquely stable and natural walking experience whether on flat, angled, smooth or uneven terrain. Kinnex 2.0 builds on this foundation with improvements in hydraulic system design, improved user experience and added product testing.



# Designed for Real Life.

- ✓ Seamlessly adjust to changes in ground slope within the first step.
- ✓ Confidently brave the elements with an IP67 certified ankle.
- ✓ Comfortably match heel heights with virtually any shoe for work or play (via App).
- ✓ Optimize energy storage and energy return with a bolt free full carbon fiber foot plate.



## Discover Kinnex 2.0

**Sophisticated microprocessor controlled resistance** that helps optimize stability whenever changes in terrain or cadence occurs.



**From sandals to barefoot in comfort**  
Kinnex accommodates heel heights from 0-2" and allows user to program up to 100 shoes.



**30° of Ankle Articulation**  
20° plantar flexion / 10° dorsi flexion allows the ankle and foot to be ideally positioned in relation to the ground.



**Smooth Rollover**  
A full length sole plate without bolts ensures seamless stance phase roll-over for superior comfort.




**Optimal Stability**  
Split keel & heel increase ground compliance and improve stability


**Water Ingress Certified**  
IP67 tested for occasional submersion in water.



**NEW Manual Lock**  
Now conveniently located inside the app or via the ankle key pad.

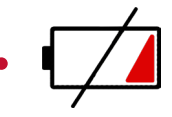


**NEW Low Battery Feature**  
Provides new enhanced notifications and improved comfort.

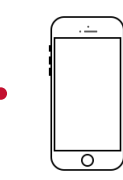

**NEW Alignment Setup Feature**  
Ankle automatically locks in a neutral position during setup for alignment needs.

**User and Prosthetist App**  
Streamlined set up, monitoring, customization and cadence reporting are a touch away with the Clinician and Patient Apps.

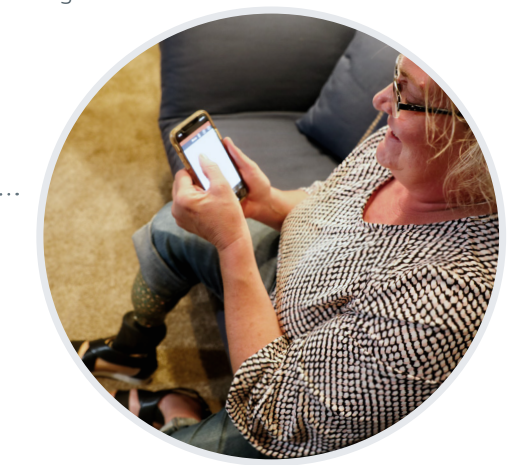

[FreedomKinnex.com](http://FreedomKinnex.com)



Discover all the ways Kinnex helps put the world at your patients' feet.