

3190T EZ ADAPT LOCK


Ultra Low Profile Modular Lock



EZ ADAPT components allow you to fabricate a low profile socket much easier, providing modular versatility to rotate as needed without any hassle.



Push Pin can be cut to various lengths

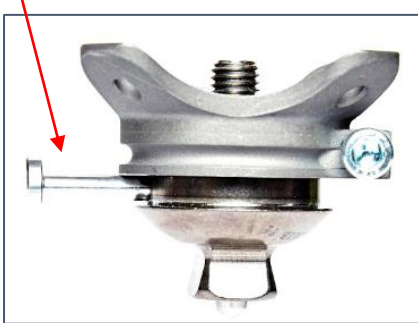

The **3190T EZLOCK** is a titanium adapter with a male pyramid on one end and threaded pylon on the opposite end. This titanium lock, with serrated pin type attachment, gives the ultimate advantage in strength and versatility when space is limited. The EZLOCK is a light space saving device which also provides the benefit of rotation. The push button can also be cut to length if needed.

The **3190T EZLOCK** is designed to work in conjunction with ST&G 3 & 4 prong rotating bases 1421S, 1221S, 1223T. Please use laminating instruction for our bases.



Fabrication Dummy

Set screws should be torqued and thread locked applied after alignment is complete.

EZLOCK Includes 2 Pin sizes (Medium & Large) and pin guide dummy

P/N	Product Description	Weight Limit	Product Weight	Build Height
3190T	Low Profile Lock	136kg / 300lbs	92g	25 mm



ST&G USA CORPORATION
2691 Saturn Street
Brea, CA 92821

TELEPHONE: (714) 524-0663
FAX: (714) 364-8113
E-MAIL: info@stngco.com

www.stngco.com